
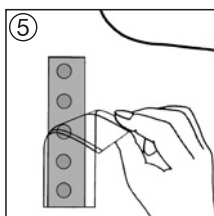
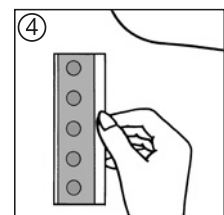
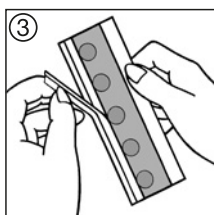
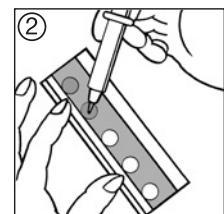
Curatest® F

Epikutantestpflaster
Emplâtre de test épicutané
Patch test strip
Apósito para pruebas epicutáneas
Adesivo para testes epicutâneos
Cerotto per test epicutanei
Epicutaaantestpleister
Epikutantestplaster
Epikutantestplâster
Epikutaanitestilaastari
Náplast pro epikutánní testy
Epikutánná testovacia náplast
Epikután próba tapasz
Plaster do testów skórnych
Obliž za epikutano testiranje
Flaster za epikutani test

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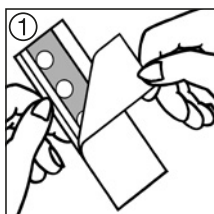en Instructions for use

Product description and features:

Transparent, waterproof strip that can be applied with test substances to perform epicutaneous tests. The strip consists of an adhesive-coated film. There are 5 test discs on the adhesive side. The strip has a backing film on its backside and two overlapping cover papers on the adhesive side.

Product composition:

Polyurethane, polyacrylate, viscose, polyethylene, polyester, silicone paper

Intended purpose:

The Curatest® F Adhesive Strip for Patch Test is intended to be used for patch-tests.

Indications:

Not applicable.

Contraindications:

Known allergy and/or hypersensitivity to any of the product components.

Side effects:

In very rare cases skin irritations and/or allergies may occur.

Application instructions:

The product is applied by healthcare professionals.

1. Preparation:

Before application, please ensure that the test area is free of hair and that the skin is dry and grease-free to ensure that the patch test strip adheres properly and that there is contact between the skin and the allergens.

2. Use:

- Remove a Curatest® F patch test strip from the folding carton and detach the wide cover paper (Fig. 1).
- Coat the test discs with the test substances (Fig. 2).
- Firmly hold the Curatest® F patch test strip by the tab and remove the folded paper strip (Fig. 3).
- Apply the Curatest® F patch test strip on the desired part of the skin with the aid of the tab. For secure fixation, gently pass over the backing film with the hand (Fig. 4).
- Peel the backing film off with the aid of the tab (Fig. 5).
- Follow-up appointment: Remove the Curatest® F patch test strip by detaching one corner of the film from the skin, and then detaching the strip by pulling steadily on the film (overstretching).

3. Disposal:

In Europe, product waste can be assigned a waste code from Chapter 18 01 and packaging waste a waste code from Chapter 15 01 of the Directive on the European Waste Catalogue (Waste Catalogue Ordinance - AVV). Packaging labelled as recyclable should be taken into the respective national recycling systems.

General instructions:

In case of a serious incident, contact the manufacturer and the responsible health authorities.

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